

Where Ideas Meet Industry



LIGHTNIN SVM Mixer Versatile Mixing for Critical Bioprocesses

LIGHTNIN's family of SanStar[™] Sanitary Vertical Mini (SVM) Mixers is designed to meet the most stringent requirements for cGMP facilities, and to provide a comprehensive array of mixer solutions for the biotech, pharmaceutical, food and cosmetics industries. These "mini" mixers have been engineered for your most critical processes, collectively providing superior performance and clean operation. Suitable for the complete range of bio-pharm and "clean" industry applications: buffer and media prep; cell culture; fermentation; product recovery; product blending and packaging.

Engineered for cleanliness in critical processes


- Minimum footprint for vessels from 20 to 500 liters.
- 300 series stainless steel seal housing and pedestal; epoxy-painted drive.
- · Designed for use in dry-running, aseptic environments.
- Vertical, in-line design for closed tank applications.
- Mixer mounts on standard 3" tri-clamp flange.
- Single dry-running seal as standard, with or without sanitary purge gland.
- Wide range of output speeds from 120 to 370 rpm, for maximum flexibility.
- Custom selections available per job for alternate speeds.

Comes with LIGHTNIN's 100% process guarantee.

- Variety of impeller configurations and variable speed drives for a wide range of process requirements; shaft length up to 45" long.
- Full clean-in-place, sterilize-in-place (CIP/SIP) features.
- Designed to comply with ASME Bioprocess Equipment Specification guidelines.





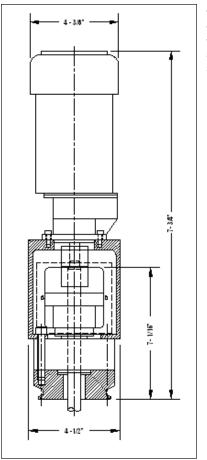

This configuration shows the standard single seal; the image on the front shows the optional double seal unit.

SVM Mixers				
Model	HP	Output Speed (RPM)	High-Efficiency Impeller Options Max. Dia., in. A510_22* A320	
SVM4S15	0.15	370	8.8"	7"
SVM7S15	0.15	245	11.2"	9"
SVM9S15	0.15	180	12.8"	10"
SVM13S15	0.15	123	12.8"	10"

 * The A510_22 is the fabricated and polished version of the A310 hydrofoil impeller.

Service

SPX Process Equipment Services

- Repair Services
- Equipment Services
- PM/PDM Services
- Additional Services
 - Installation and Startup
 - Maintenance and Repair
 - Asset Management

For 24-hour response, call 1-888-MIX-BEST (U.S. and Canada)

Call the SPX Process Equipment Experts

For more information about the LIGHTNIN SVM Mixers or to discuss your application, contact your local SPX Process Equipment Sales Representative. Or visit our website at www.spxpe.com.

Or call 1-888-MIX-BEST (U.S. and Canada)



PROCESS EQUIPMENT



SPX Process Equipment, LIGHTNIN Operation 135 Mt. Read Blvd., Rochester, NY 14611 Phone: (888) 649-2378 (MIX-BEST), (U.S. and Canada) or +1 (585) 436-5550 E-mail: LIGHTNIN@spx.com

For more information about our worldwide locations, approvals, certifications, and local representatives, please visit www.spxpe.com.

SPX Corporation reserves the right to incorporate our latest design and material changes without notice or obligation. Design features, materials of construction and dimensional data, as described in this bulletin, are provided for your information only and should not be relied upon unless confirmed in writing.

Issued: 09/2008 B-914

Copyright © 2008 SPX Corporation

Your local contact:

when you purchase a LIGHTNIN SVM Mixer.

QAD for cGMP Requirements

To demonstrate that FDA-mandated Good Manufacturing

Practices (cGMP) guidelines have been met, you receive

full Quality Assurance Documentation (QAD) as standard

as Standard



Global Headquarters 13515 Ballantyne Corporate Place Charlotte, North Carolina 28277 United States Where Ideas Meet Industry

Where Ideas Meet Industry

www.spx.com